
# Kit Cãnula Pulse Block

Reg. ANVISA: 80195529005.
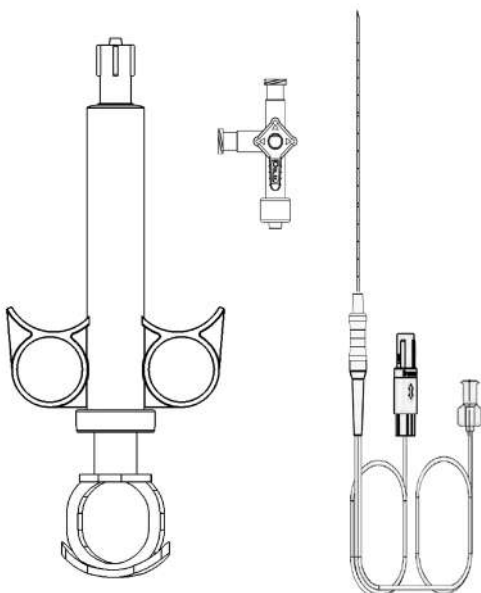
## BRAMSYS

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AV. DR. ALBERTO SAN JUAN, 350

18530-000 – TIETÊ – SP

CNPJ: 04.729.995/0001-87

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## KIT CÂNULA PULSE BLOCK

### FINALIDADE:

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### **PULSE BLOCK CANNULA KIT**

#### **PURPOSE:**

The Pulse Block Cannula Kit are intended for performing plexus and peripheral nerve block procedures.

#### **INDICATIONS OF USE**

The Pulse Block Cannula Kit is indicated for Plexus and Peripheral Nerve Block procedures, as well as general or conductive anesthesia for neurolytic block, cranial nerve block, celiac plexus anesthetic block, nerve plexus anesthetic block (lumbosacral, brachial , cervical) for pain management, Lumbar sympathetic anesthetic block, Sympathetic anesthetic block, Temporomandibular joint block, Local anesthetic stellate ganglion block, Neurolytic stellate ganglion block, Peripheral nerve block, Peripheral nerve block - Anesthetic blocks of nerves and neurovascular stimuli, Block of the autonomic nervous system, Paraspinous facet block, Phenolic, Alcoholic or Botulinum Toxin Block - By Body Segment, Neurolytic block of cranial or cervicothoracic nerves, Neurolytic block of the celiac plexus, lumbar or thoracic sympathetic nerve, Neuroma extirpation, Intraneural or intrafascicular microneurolysis of two or more nerves, Intraneural or intrafascicular microneurolysis of a nerve, Multiple microneurolysis, Single microneurolysis, Presacral or genitofemoral nerve neurectomy, Neurolysis of compression syndromes, Diagnostic or therapeutic joint puncture (infiltration) - Guided or not by imaging method, Diagnostic or therapeutic extra-articular puncture (infiltration/dry needling) - Guided or not by imaging method, Radioscopy for follow-up of surgical procedures and Microsurgical treatment of compressive neuropathies, with the aid of percutaneous electrical stimulation for target location. These procedures can be performed together, or separately, at the discretion of the physician.

## COMPOSITION

(1) blocking cannula, (1) graduated control injector, (1) 3-way adapter, (1) surgical planning scale, (1) surgical marker.

## USE INSTRUCTIONS

- Make sure the packaging is not damaged.
- Remove Blocking Cannula, Graduated Control Injector, 3-Way Adapter, Surgical Planning Scale and Surgical Marker from sterile packaging in a clean environment using aseptic techniques.
- Make sure that all items: Blocking Cannula, Graduated Control Injector, 3-Way Adapter, Surgical Planning Scale and Surgical Marker are intact.
- Position the patient according to the surgical procedure to be performed.

NOTE: If necessary, use the Surgical Planning Scale and the Surgical Marker to assist in puncture planning and the desired target location.

- Prepare skin puncture area using aseptic and antiseptic techniques
- Perform the percutaneous puncture of the Blocking Cannula in the patient.

NOTE: If necessary, use fluoroscopy, ultrasound, or x-ray to aid puncture and target monitoring.

- Connect the stimulation/neuromodulation cable to the BMS-50N Generator Unit to help locate the nerve through the stimulus and, if necessary, perform temperature-controlled neuromodulation.

NOTE: If necessary, use 3 way Adapter to collect and prepare the substance that will be administered during the procedure.

- Connect the graduated control injector containing the desired solution to the injector tube with luer lock connection to perform the procedure.
- After the end of the procedure, all materials must be de-characterized and discarded in an appropriate place, following the regulations applicable to the disposal of sharps and contaminated material adopted by the hospital.

## CONTRAINDICATIONS

- Infection or hematoma near the site chosen for the puncture;
- Neurological alterations known in the preoperative period are a relative contraindication;
- Stimulation: do not perform during defibrillation or cardioversion;
- Changes in coagulation, eg during anticoagulant therapy are a relative contraindication;
- Lymphangitis (for example, in case of axillary access);
- Post-Lymphadenectomy status (eg, in case of axillary access);
- Acute plexus neuritis;
- Inability to obtain patient consent;
- Known hypersensitivity to the materials used.

## STORAGE CONDITIONS

The Pulse Block Cannula Kit must always be stored protected against all types of moisture, dirt and contact with flammable, corrosive, toxic or explosive substances. The packaging must not be punctured or torn, it must be kept away from sharp or pointed instruments.

Do not place it under weight or volumes, as they may damage its technical characteristics and compromise its integrity.

Do not store the product where it can be exposed to direct sunlight, high temperatures, humidity or dust as this will compromise its integrity/quality.

Do not place them next to or together with different categories, as they may damage their technical characteristics and compromise their integrity / quality.

Storage temperature: 5°C to 60°C

Storage humidity: 10% to 80%

## PRECAUTIONS, RESTRICTIONS AND WARNINGS

- The user, before using the product, must carefully read the instructions for use that will help correct and safe use, as well as protecting against possible risks.
- Do not use the product after its expiration date, always check if the sterilization indication is correct.
- The product must only be used by qualified medical professionals, who have the knowledge and experience necessary for the proper use of the product, knowledge of handling sterile products and, mainly, of the importance of maintaining a totally sterile environment in the operating rooms. , where care and precaution with every detail are essential so that there is no risk of contamination.

- Check the integrity of the packaging before using the product, in case of violation, DO NOT USE, as the guarantee of its sterility will be compromised.
- This product is sterilized by Ethylene Oxide (ETO).
- The product must be transported and stored, always protected against all types of moisture, dirt and contact with flammable, corrosive, toxic or explosive substances.
- The packaging must not be punctured or torn, it must be kept away from sharp or pointed instruments.
- Do not store or store this product next to products of different categories.
- It is prohibited to reprocess this product.
- The product must not be reused and/or re-sterilized, it is for SINGLE use, and must be destroyed and disposed of in HOSPITAL WASTE.

This product is totally safe if used by qualified professional and respecting the usage specifications described in this document.

**MS Registration:** 80195529005  
**Technical Responsible:** Marcos Rojas  
**CREA:** 5061214879

### **WARRANTY CERTIFICATE**

**BRAMSYS INDÚSTRIA E COMÉRCIO LTDA**, properly authorized by the National Agency of Sanitary Surveillance under number UY44WX94W79M, guarantees for **24 MONTHS**, starting from the sterilization date, the quality and the sterile condition of this product, except in the cases of bad usage or noncompliance of the orientations written in this document.

### **CUSTOMER SERVICE (CS)**

Phone: +55 11 5013 0010  
E-mail: [sac@bramsys.com.br](mailto:sac@bramsys.com.br)  
Site: [www.bramsys.com.br](http://www.bramsys.com.br)

### **MODELOS / MODELS**

KPB1-1805-5; KPB2-1805-5; KPB1-1810-5; KPB1-1810-10; KPB2-1810-5; KPB2-1810-10;  
KPB1-1815-5; KPB1-1815-10; KPB2-1815-5; KPB2-1815-10; KPB1-2005-5; KPB2-2005-5;  
KPB1-2010-5; KPB1-2010-10; KPB2-2010-5; KPB2-2010-10; KPB1-2015-5; KPB1-2015-10;  
KPB2-2015-5; KPB2-2015-10; KPB1-2105-5; KPB2-2105-5; KPB1-2110-5; KPB1-2110-10;  
KPB2-2110-5; KPB2-2110-10; KPB1-2115-5; KPB1-2115-10; KPB2-2115-5; KPB2-2115-10.